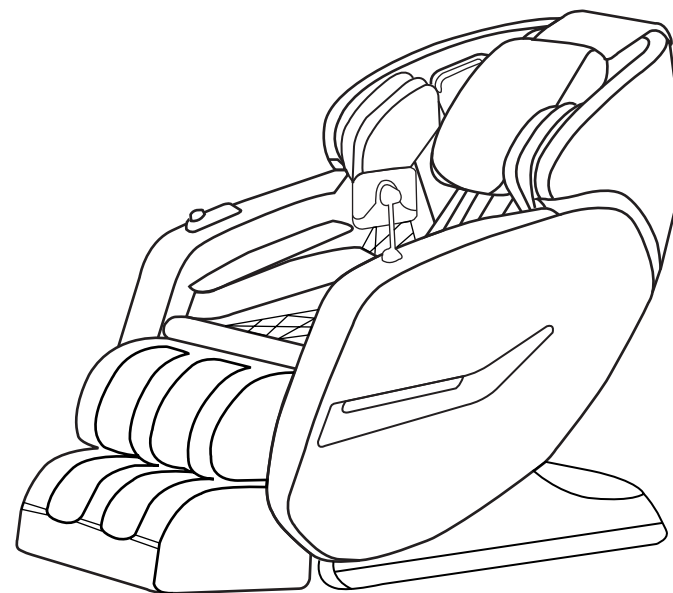


# **Luxury multifunctional massage chair**



**User's manual**

Thank you for choosing this product. Before you use it, please read this instruction manual in order to master the correct use of this product. When referring to it, please pay special attention to important safety warnings. At the same time, keep the instruction manual for future reference.  
Note: The company reserves the right to change and explain the design of the product. If there is any change without notice, the color of the product is subject to the actual product.

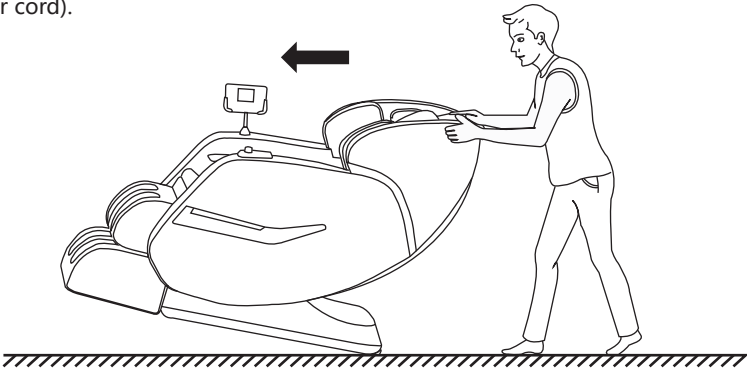
## Contents

Safety and maintenance matters.....	1-2
Use safety warning .....	3
Part Name .....	4
Product Features .....	5
Hand control introduction.....	6
Button function of hand control.....	7
Instructions .....	8-9
Product specifications.....	9

## Instructions

### 3. Moving method

After confirming that all the wires are pulled out of the massage chair, press the backrest down to a certain angle (when the center of gravity falls on the pulley), you can push the massage chair forward with your hand, and then move the massage chair to a suitable position. Back to normal angle. (Note: When moving, be sure to turn off the power switch before unplugging the power cord).



#### warn

The moving pulley behind the massage chair is only suitable for short-distance and slow-speed pushing. For long-distance transportation, please use other auxiliary moving tools to move the massage chair!

## Product specifications

- Product name: luxury multifunctional massage chair
- Rated voltage: 220V~
- Rated frequency: 50Hz
- Rated power: 120W
- Rated working time: 15-30min
- Executive standard: GB4706.10-2008/ GB4706.1-2005
- Product safety design: Class I

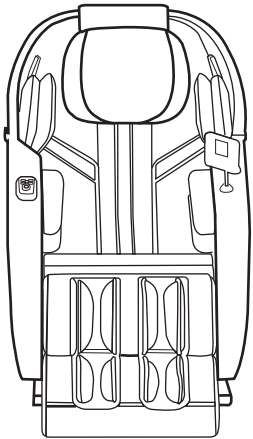
## Instructions

### 1. Bluetooth audio playback operation instructions


- (1) Bluetooth is turned on by default when it is turned on and voice prompts;
- (2) Turn on the Bluetooth function of the mobile phone, find the Bluetooth connection of this machine, link password 8888, pairing is successful, turn on the music and start playing.

### 2. Handrail disassembly



When the doorway is too small and the chair cannot be moved into the door, the armrest can be removed and transported separately.



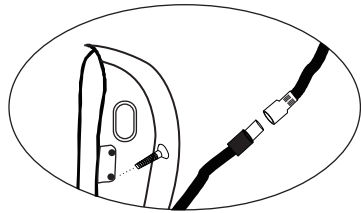1. Unscrew the hexagon socket screws on both sides
2. Unplug the speaker connector



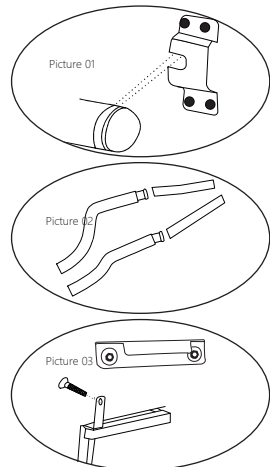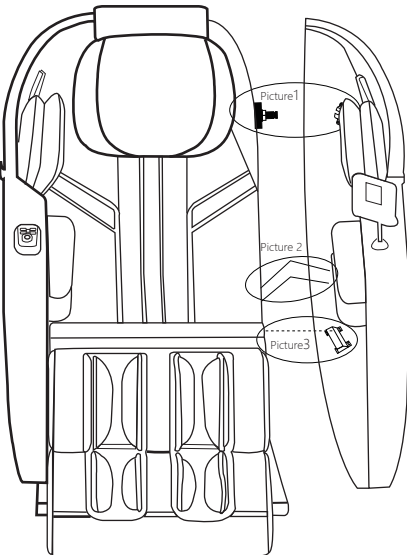1. Connect the trachea in Figure 2 according to the color;
2. Then fit the two card slots of Figure 1 and Figure 3 together. When the two card slots are installed, the air should be managed smoothly and bend carefully;
3. Finally, lock the screws in Figure 3;



## Safety and maintenance matters

### 1. Important safety warning

- Do not allow children to approach the movable part of this appliance;
- Please use a well-grounded power source that matches this appliance;
- Please unplug the power plug after use or before cleaning, so as not to cause equipment or personal damage;
- Please use this appliance according to the instructions in this manual;
- Do not use unrecommended accessories and accessories;
- Do not use this appliance outdoors;
- Please read this operation manual carefully before starting to use;
- It is forbidden to use this appliance for other purposes not listed in this manual;
- The recommended use time is 20 minutes;
- Do not use this appliance when the cover cloth is broken or damaged;
- Do not use this appliance when the outer shell or outer cover of this appliance is damaged;
- Do not use any objects to cover the cooling holes of this appliance;
- Do not let any objects fall into this appliance;
- Do not sleep during the use of this appliance;
- Do not use this appliance when you are drunk or feel unwell;
- Do not use this appliance within one hour after eating;
- Do not use excessive force when using this appliance to avoid injury;
- It is forbidden to carry fragile, explosive and flammable items such as mobile phones and lighters in the massage and compressed areas.

### 2. Use environment

- Do not use this appliance in a highly humid environment;
- Do not use this appliance immediately when the ambient temperature changes sharply;
- Do not use this appliance in an environment with severe dust or corrosive gas;
- Do not use this appliance in a small space or a place with poor air circulation;
- Normal room temperature 0C to 35°C;
- At a relative humidity of 20% to 80%, do not use this appliance in an excessively humid room (such as a bathroom);
- Do not use this appliance outdoors;
- Avoid direct sunlight;
- Keep this appliance away from heat sources (such as the side of the stove, the air outlet of the air conditioner);
- Do not use this appliance in thunder and lightning;
- Do not use this appliance in places with excessive dust.

### 3. Those who are not suitable to use massage chairs

- Patients with osteoporosis;
- Those who have heart problems and wear electronic medical devices such as pacemakers;
- Those who have a high fever;
- Pregnant women or menstruating women;
- People with physical trauma or body surface disease;
- Children under the age of 14 and unconscious persons are prohibited from using this product when unsupervised;
- Those who need to recuperate or feel unwell after being ordered by a doctor;

## Safety and maintenance matters

- Do not use this appliance when the body is wet;
- The appliance has a heating surface, so people who are not sensitive to heat must be careful when using it.

### 4. Safety matters

- Check whether the voltage meets the specifications of this appliance;
- Do not use wet hands to pull or plug the power plug;
- Do not allow water to enter this appliance to avoid electric shock or damage to this appliance;
- Do not pull the power cord when plugging or unplugging the power plug, and don't be too rough;
- Do not damage the wire or tamper with the circuit of this appliance;
- Do not use a damp cloth to clean the live parts such as the power switch and plug of this appliance;
- In the event of a power failure, you should leave the appliance to avoid injury when the power is suddenly restored;
- If you feel that this appliance is abnormal during use, please stop using it immediately and consult your local supplier;
- If you feel an abnormality during the use of this device, please stop using it immediately and consult a healthcare practitioner;
- This device is not intended to be used by people (including children) who are physically disabled, sensory or neurologically defective, or lack experience and common sense, unless they are supervised and guided by those responsible for their safety on how to use the device. Children must be supervised to ensure that they do not use equipment to play with;
- If the power cord is damaged, in order to avoid danger, it must be replaced by professionals from the manufacturer, its maintenance department or similar departments.

### 5. Product care and maintenance

- This appliance can only be repaired by the service center designated by the manufacturer, and the user cannot disassemble or repair it by himself;
- Please remember to turn off the main power switch after use;
- Do not use this appliance when the power socket is loose;
- If the appliance will not be used for a long time, the power cord of this appliance should be rolled up and stored in a dry and dust-free environment;
- Do not store this appliance near high temperature or open flame, and do not expose it to direct sunlight for a long time;
- Please clean this appliance with a thousand cloth, do not use thinner, benzene or alcohol;
- The mechanical device of this appliance is specially designed and manufactured without special maintenance;
- Do not pierce this device with sharp objects;
- Do not roll or drag this appliance on uneven ground, it needs to be lifted and moved;
- Please use it intermittently. Do not allow this appliance to run continuously for a long time.

### 6. Treatment of common phenomena

- Motor sound during operation: it is a normal operation sound;
- The controller cannot work normally: Check whether the power plug and the power socket are connected reliably; check whether the power switch is turned on;
- The operation of the appliance stops: when the rated time is up, the appliance will automatically activate the shutdown button; if the appliance has been working continuously for too long, the temperature control protection switch will automatically shut down, allowing the appliance to rest for half an hour before using it.

### 7. Storage conditions

- Storage room temperature -5°C~35°C;
- Relative humidity 20% to 80%, no corrosive gas and well ventilated indoor.

## Button function of hand control











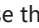
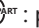

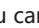






### Boot up

- Standby key: Press the "⏻" key to turn on/off the massage function of the whole machine.
- When the power is turned on, press this button to lie down to the initial angle. After completion, the first automatic program "President Health Care" will run by default. The default working time is 20 minutes. Press this button again and wait for another 20 seconds. Turn off all functions of the massage chair.

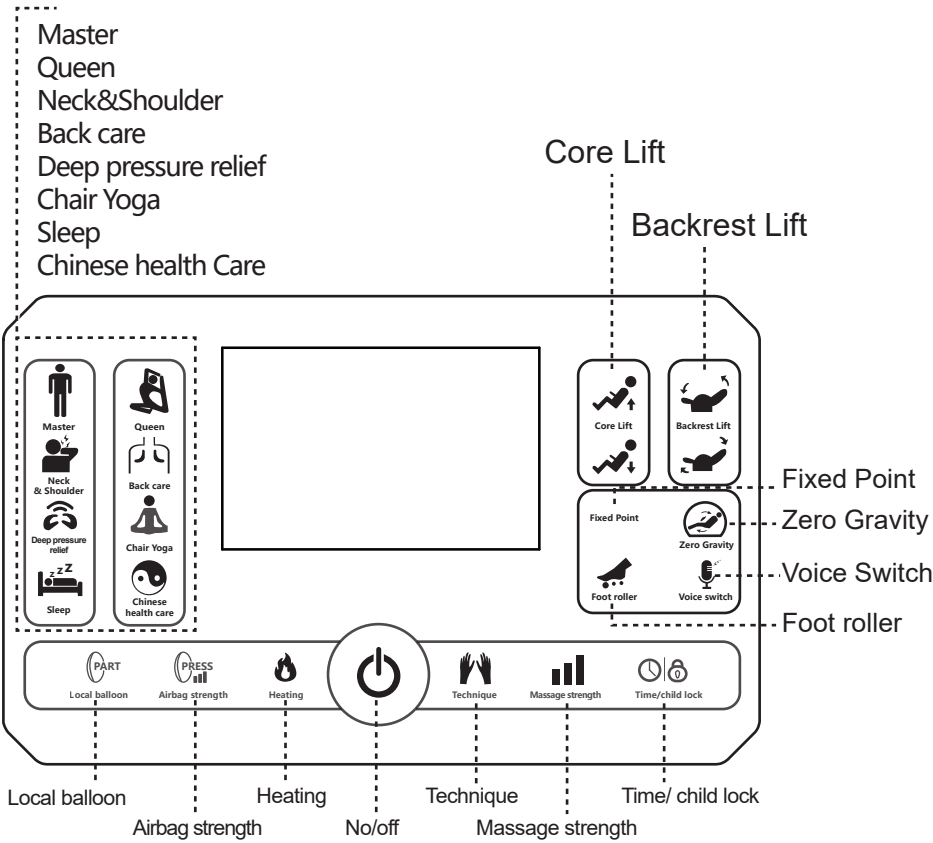
### Automatic mode

- Automatic massage program key: According to your body condition and massage time, you can choose: Master \Queen\Neck&Shoulder\Back care\Deep pressure relief\Chair Yoga Sleep \Chinese health Care, any of the eight massage programs for targeted massage. (See LCD screen instructions for more information), you can also click the shortcut button AUTO key to select in turn.


### Other functions

- Manipulation : press the "Manipulation" key to select 6 massage techniques: kneading, tapping, kneading, massage, tapping and finger pressing.
- Movement up and down: long press the "" key to adjust the massage of the upward part of the manipulator and long press the "" key to adjust the massage of the downward part of the manipulator. (shortcut keys "انر" and "انر" can also be used to operate movement lifting and lowering)
- Heating: press the "" key to turn on/off the waist heating function. (You can also use the shortcut key")
- Foot roller: press the "" key to turn on/off the foot roller;
- Time : press the "time" key to select 10 minutes /20 minutes /30 minutes for regular massage;
- child lock : press the "child lock - zero gravity : press the "zero gravity" key to turn on/off zero gravity, press the "" key for a long time to raise the backrest, and press the "" key for a long time to lower the backrest; (You can also use shortcut keys)
- Airbag part : press the "Airbag part" key to massage the arms, shoulders, calves, etc. of the part of the air bag;
- Airbag strength : Press the "Airbag Strength" key to adjust the massage strength of the air bag by 1-3 steps; (You can also press the shortcut button "" to adjust the force of the airbag in 1-3rd gear);
- Fixed-point/local massage: Press the "Fixed-point/local" key to adjust the massage parts of the manipulator, such as the interval massage of neck, shoulder, back, waist, buttocks, thighs, etc.;
- Shoulder height adjustment: When the shoulder height is not in place during massage, you can double-click the "" key to adjust the upward movement of the manipulator, and double-click the "" key to adjust the downward movement of the manipulator until it is suitable for your shoulder height;
- Voice switch : press the "Voice switch - Massage intensity : press the "Massage intensity 

# Hand control introduction




# Shortcut button

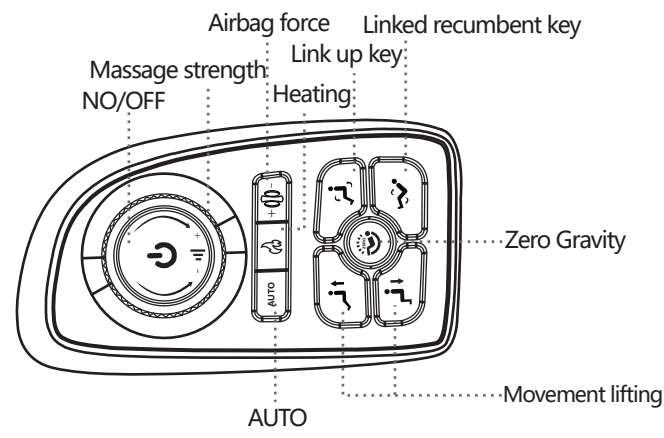

# Use safety warning


## ! warn

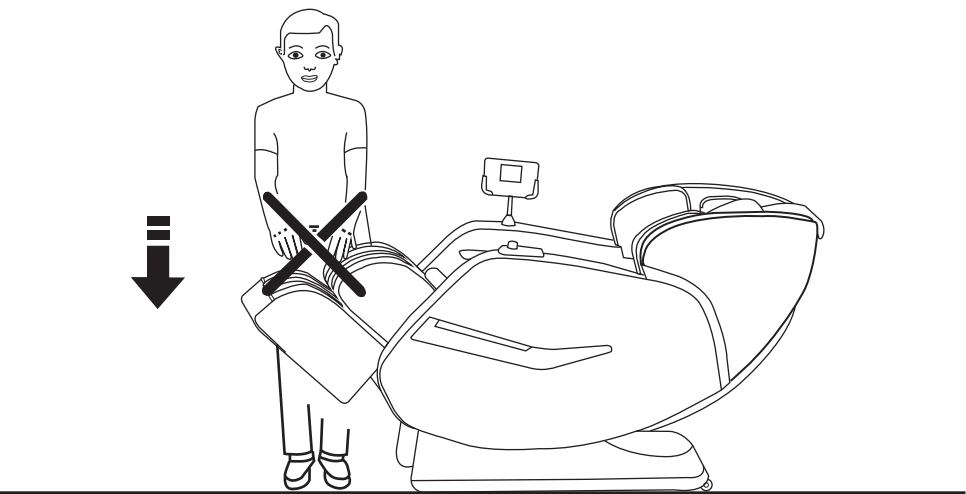When the calf frame of the massage chair rises, people must not apply gravity to the calf frame to prevent the massage chair from falling forward to avoid product damage or personal injury.




## ! warn

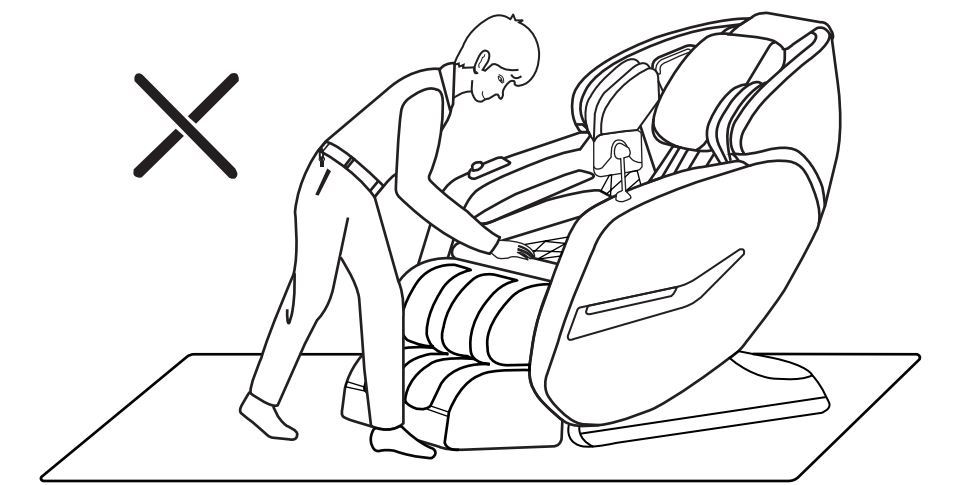Before sitting on the massage chair, it is necessary to confirm that the massage manipulator cannot stay in the cushion position, so as to avoid discomfort to the buttocks and damage to the massage manipulator during the downward pressure of the human body.



## Part Name

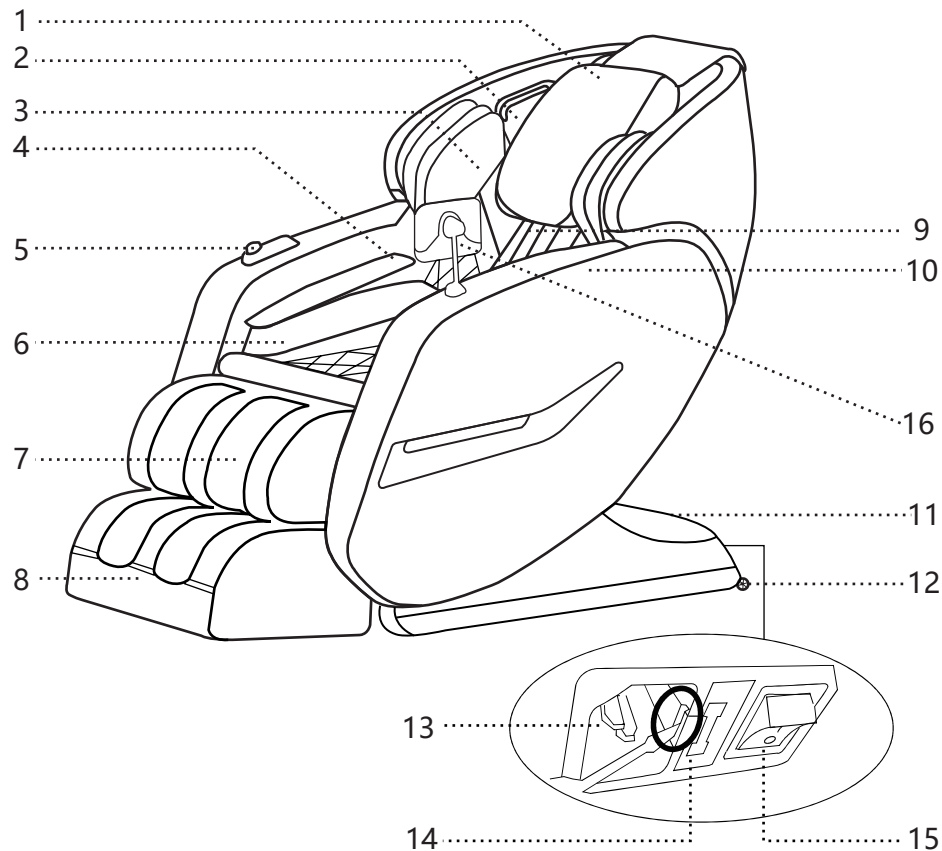

- |                             |                          |                                 |
|-----------------------------|--------------------------|---------------------------------|
| 1, Head cushion             | 2, 3D digital audio      | 3, Shoulder pneumatic component |
| 4, Arm pneumatic component  | 5. Shortcut Button Keys  | 6, Back cushion                 |
| 7, Calf Pneumatic Component | 8, Foot Roller           | 9. Armrest                      |
| 10, Armrest                 | 11, Power box back cover |                                 |
| 12, Moving wheel            | 13. Power cord jack      | 14, Fuse                        |
| 15, Power switch            |                          |                                 |
| 16, Large screen controller |                          |                                 |

## Product Features

- This product is designed with a set of 3D manipulators that can walk up and down, stretch back and forth, and mute design.
- Designed with automatic shoulder position detection and micro adjustment function; The automatic detection function of human body curve and massage finger pressure points can automatically adjust the front and rear extension of the manipulator according to the detected human body curve and massage finger pressure points, making the massage more humanized and scientific.
- Featured automatic massage functions: Master, Queen, Neck& Shoulder, Back care, Deep pressure relief, Chair Yoga, Sleep, and Chinesehealth care.
- Manual massage range selection of manipulator: including fixed point, upper back, lower back, local and whole massage range selection.
- It is equipped with manual massage function, including six massage methods: stop, knead, pat, knead and pat, finger press, and knock.
- It is equipped with massage adjustment functions, including 6 levels of massage intensity adjustment, 3 levels of air pressure intensity adjustment and other massage adjustment functions.
- It has upper body air pressure massage function: palm air pressure massage function, shoulder air pressure massage function, arm air pressure massage function, and three levels of air pressure intensity adjustment.
- It has lower body air pressure massage function: lower leg air pressure massage function, foot air pressure massage function, and three levels of air pressure intensity adjustment.
- The guide rail of the backrest frame adopts the curve guide rail technology, which fully conforms to the curve of the back of the human body, and can enjoy the intimate massage from the head to the thigh.
- Back heating function: far infrared carbon fiber heat source is used.
- Music playing system with built-in 3D digital audio on the left and right. The Bluetooth function of the intelligent device can be used to play music wirelessly.
- It is equipped with backrest linkage lifting adjustment function and zero gravity automatic reclining angle adjustment function.
- It is equipped with foot roller massage function, and the speed is adjustable in three gears.
- The lower leg is equipped with elastic yellow expansion and adjustment function. The extended length of the foot is determined according to the length of the user's leg. The foot is retracted by the spring force.
- Automatic sit up function after shutdown, and automatic reset function of manipulator.
- Large screen manual controller and shortcut button are configured for more convenient operation.